


# Metal Chill Renk Deęiřtiren Soęutucu Rack

- 4 °C altında tahmini süre: 2,5 saat.
- Kullanımdan önce -20 derecede 8 saat dondurulur.
- Isıyla renk deęiřtirebilen organik boya ile boyanmıřtır.
- Isı yükseldikçe koyudan aıęa renk deęiřir.
- Yüksek kalitede alüminyum materyalden üretilmiřtir.
- Sterilizasyon sıvıları ile temizlenebilir.

## > NB- Metal Chill01

- 40 x 1.5-2.0ml tüp
- 4 x 0.5 ml tüp
- 96 x 0.2ml Plate

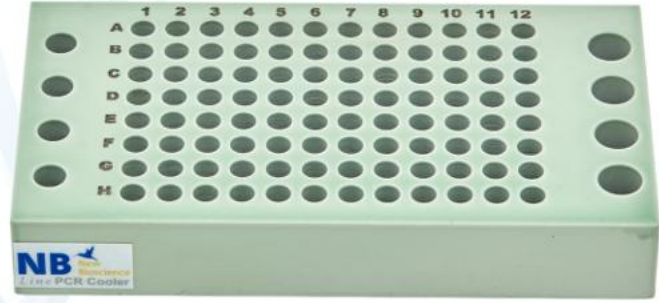

## > NB- Metal Chill02

- 4 x 1.5-2.0ml tüp
- 4 x 0.5 ml tüp
- 96 x 0.1 ml Plate 4'lü strip tüplere uygun



+90 312 220 22 80

bm@bmlabosis.com

Kızılırmak Mah. 1441. Cad. Meva İş Merkezi No:41/8 A


www.bmlabosis.com

# Metal Chill Color Changing Cooling Rack

- Estimated time under 4 °C: 2.5 hours.
- Must be frozen at -20 °C for 8 hours before use.
- Coated with heat-sensitive organic dye.
- Color changes from dark to light as temperature increases.
- Made from high-quality aluminum material.
- Can be cleaned with sterilization liquids.


## > NB- Metal Chill01

- 4 x 1.5-2.0 ml tube
- 4 x 0.5 ml tube
- 96 x 0.2 ml plate

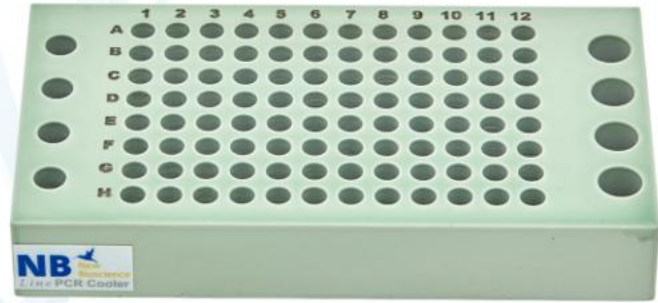

## > NB- Metal Chill02

- 4 x 1.5-2.0 ml tube
- 4 x 0.5 ml tube
- 96 x 0.1 ml plate suitable for 4-strip tubes

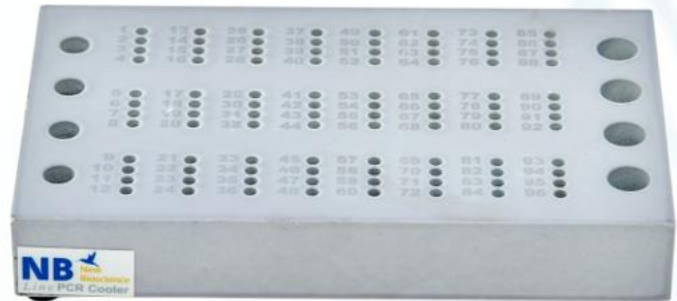

+90 312 220 22 80

bm@bmlabosis.com

Kızılırmak Mah. 1441. Cad. Meva İş Merkezi No:41/8 A

www.bmlabosis.com