

## AUTOMATIC NUCLEIC ACID PURIFICATION SYSTEM





- +90 312 220 22 80
- mww.bmlabosis.com
- ≥ bm@bmlabosis.com
- ♥ Kızılırmak Mah. 1441. Cad. Meva İş Merkezi No:41/8 Ankara/TÜRKİYE







## Product Features



Industrialized control system makes stable operation for 24 hours

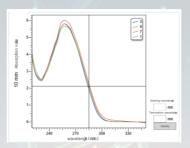

High product yield and good purity

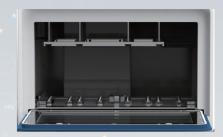

Automatic nucleic acid extraction processing can be performed on 32/96 samples simultaneously, greatly liberating the hands of researchers



Supporting reagents can be applied to various samples such as swabs, serum plasma, tissues, plants, whole blood, fecal multiple specifications of single/16T/32T/48T/96T soil, bacteria, etc., and have to meet different customer needs




Self developed intelligent operation software and touch screen makes operation convenient and fast




Disposable sheath insulates magnetic rods and samples, and machine is equipped with UV sterilization and air filtration systems to reject cross contamination adsorption




## Application Case

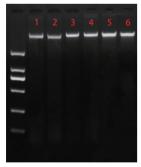

Electrophoresis test results of faeces and soil samples after extraction 1,2: human faeces 3,4: cat faeces 5,6: soil

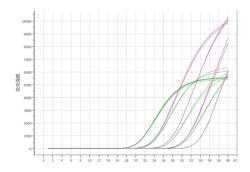

UU sample extracted qPCR analysis results

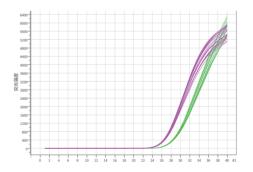

NG sample extracted qPCR analysis results (including internal standard)

No.	type	potency	unit	A260	A280	260/280	260/230	sample
1	RNA	556.505	μg/ml	13.913	6.636	2.097	2.393	- spleen
2	RNA	540.713	μg/ml	13.518	6.441	2.099	2.079	
3	RNA	799.469	μg/ml	19.987	9.558	2.091	2.352	kidney
4	RNA	847.294	μg/ml	21.182	10.133	2.090	2.269	
5	RNA	1087.187	μg/ml	27.180	12.870	2.112	2.344	liver
6	RNA	980.632	μg/ml	24.516	11.626	2.109	2.329	liver



## > Recommend Kits

PRODUCT NAME	PACKING	CAT. NO.	
Magpure animal tissue genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP01R/BFMP01R96	
Magpure Zebra fin tissue genomic DNA purification kit (pre-filled package)	32Т/96Т	BFMP01RD/BFMP01RD96	
Magpure whole blood genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP02R/BFMP02R96	
Magpure plant genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP03R/BFMP03R96	
Magpure plant genomic DNA purification kit (polysaccharide and polyphenol)	32Т/96Т	BFMP03RP/BFMP03RP96	
Magpure virus DNA purification kit(pre-filled package)	32Т/96Т	BFMP04R/BFMP04R96	
Magpure dried blood spots genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP05R/BFMP05R96	
Magpure oral swab genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP06R/BFMP06R96	
Magpure total RNA purification kit(pre-filled package)	32Т/96Т	BFMP07R/BFMP07R96	
Magpure virus DNA/RNA purification kit(pre-filled package)	32Т/96Т	BFMP08R/BFMP08R96	
Magpure plasmid DNA purification kit(pre-filled package)	32Т/96Т	BFMP09R/BFMP09R96	
Magpure bacteria genomic DNA purification kit (pre-filled package)	32Т/96Т	BFMP10R/BFMP10R96	
Magpure Plasma DNA purification kit (pre-filled package)	32Т/96Т	BFMP11R/BFMP11R96	
Magpure FFPE genomic DNA purification kit(pre-filled package)	32Т/96Т	BFMP12R/BFMP12R96	
Magpure virus RNA purification kit(pre-filled package)	32Т/96Т	BFMP14R/BFMP14R96	
Magpure soil and faeces genomic DNA purification kit (pre-filled package)	32Т/96Т	BFMP15R/BFMP15R96	
96 deep well plate (2.2ml) (U-shaped)	96/	BFMH01	
8-strip tips	100/	BFMH02	
96 deep well plate (2.2ml) (V-shaped)	50/	BFMH07	
96-strip tips	50/	ВҒМНО8Е	







MODEL	BFEX-32E	BFEX-96E		
Reagent type	Magnetic beads method			
Processing volume	20~1000 L			
Sample throughput	1~32	1~96		
Magnetic rod origin calibration	Software calibration (without dismantling)			
Operation screen	10.1 inch LCD touch			
Operation system	7×24 hours industrial			
Program management	New, edit and delete programs,input and output programs freely			
Uniformity of wells	CV<3%			
Heating temp. range	RT.+5~125			
Working temp.	10~40			
Working humidity	≤90%			
Mixing speed	Multi adjustable			
Multi absorption modes	Normal/Step/Cycle			
Retention of magnetic beads	≥99%			
Disinfection	UV lamp			
Interfaces	USB			
Power supply	AC110~240V 50/60Hz			
Size (mm)	340×360×395	616×375×394		
Net Weight (kg)	24	30		
Filter system	HEPA filter system(Exhaust type)			
Built-in lighting	Standart			
Testing of tips	N/A Standart			
Code scanner	Optional			
Antidripping liquid	N/A	Optional		